Photron



Compact high-speed camera system

For use with a wide range of general scientific and industrial applications the Photron FASTCAM Mini UX high-speed camera provides outstanding imaging performance at a very attractive price performance ratio.

Two Mini UX camera models provide 1.3-megapixel (1280 x 1024 pixels) image resolution with frame rates up to 2,000fps from the Mini UX50 and 4,000fps from the Mini UX100. Both models are available with recording memory options up to 32GB providing extended recording times and triggering flexibility.

Using innovative proprietary CMOS image sensor technology, the FASTCAM Mini UX achieves high light sensitivity from a small image sensor (10µm pixel pitch) through the utilization of microlenses to increase effective Fill Factor. At maximum image resolution the image sensor is fully compatible with readily available 1-inch C-mount lenses offering a wide choice of small, light weight, rugged and high aperture objective lenses.

The FASTCAM Mini UX features a rugged design suitable for operation in high shock and vibration environments and a compact camera body (120mm x 120mm x 93mm) weighing just 1.5kg. This small and rugged camera design makes the FASTCAM Mini UX ideally suited to on-board and off-board automotive safety testing and many other applications where a compact size and compatibility with standard optical systems is required.

Standard operational features of the FASTCAM Mini UX include a Gigabit Ethernet Interface for reliable system control with high-speed data transfer to PC, and the ability to remotely switch off cooling fans to eliminate vibrations when recording at high magnifications.

FASTCAM MiniUX

Model UX50 / UX100

1.3-Megapixel CMOS Sensor: Mini UX50

1280 x 1024 pixels at 2,000fps 1280 x 800 pixels at 2,500fps **Mini UX100**

1280 x 1024 pixels at 4,000fps 1280 x 1000 pixels at 5,000fps 1280 x 800 pixles at 6,250fps

Maximum Frame Rate:

160,000fps (Mini UX50 type 160K) 204,800fps (Mini UX100 type 200K) 800,000fps (Mini UX100 type 800K)

Class Leading Light Sensitivity: ISO 12232 Ssat

ISO 10,000 monochrome

• ISO 5,000 color

Global Electronic Shutter:

Minimum Shutter speed 3.9µs (to 1µs dependent on frame rate selection)

Dynamic Range (ADC):

12-bit monochrome, 36-bit color

Compact and Lightweight:

120mm (H) x 120mm (W) x 93mm (D) 4.72" (H) x 4.72" (W) x 3.66" (D) Weight: 1.5Kg (3.13 lbs.)

Internal Recording Memory: 4GB, 8GB, 16GB, or 32GB

Fast Gigabit Ethernet Interface:

Provides high-speed image download to a standard notebook/PC

Flexible Frame Synchronization:

Frame rate may be synchronized to external unstable frequencies

1-Inch C-mount Compatible Sensor Size:

Also supplied with integrated Nikon G-type lens mounts.

High-G Rated:

Suitable for application in high-G environments; operation tested to 100G, 10ms, 6-axes



Light Sensitivity:

Expressions of light sensitivity in high-speed cameras can be confusing as a variety of differing measurement techniques are used. Photron publishes light sensitivity figures for its products using the ISO 12232 Ssat Standard.

FASTCAM MINI UX	ISO 12232 Ssat
Monochrome models	ISO 10,000
Color models	ISO 5,000

ISO 12232 Ssat values published by Photron for both monochrome and color cameras are measured excluding infrared sensitivity as defined by the ISO standard measurement procedure ISO 14524.

Monochrome sensors used in the FASTCAM Mini UX cameras are supplied without an IR absorbing filter, extending the camera spectral response beyond 900nm. When the sensitivity of the FASTCAM Mini UX camera is measured to tungsten light including near IR response an equivalent value of ISO 25,000 is obtained.


Image Sensor:

The FASTCAM Mini UX system uses an advanced CMOS image sensor optimized for light sensitivity and high image quality that is unique to Photron.

A 10-micron pixel pitch gives a sensor size at full image resolution of 12.8 x 10.24mm (diagonal 16.39mm).

Lenses designed for both 1-inch C-mount and FX / DX (APS-C digital SLR) formats are compatible with the FASTCAM Mini UX at full image resolution.

Sensor Type	Proprietary Design Advanced CMOS				
Maximum Reolution (pixels)	1280 x 1024 pixels				
Sensor Size / Diagonal	12.80 (H) x 10.24mm (V) / 16.39mm (D)				
Pixel Size (microns)	10µm x 10µm				
Quantum Efficiency	62.6% at 630nm				
Fill Factor	80%				
Color Matrix	Bayer CFA (single sensor)				
ISO 12232 Ssat sensitivity	ISO 10,000 monochrome ISO 5,000 color (monochrome sensor equivalent ISO 25,000 including near IR response)				
Shutter	Global Electronic Shutter 3.9µs up to 1.01µs at maximum frame rates				

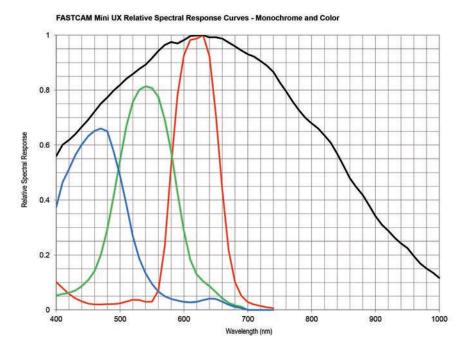

Model	Mini UX50	Mini UX100				
Full Frame Performance	2,000fps	4,000fps				
Maximum Frame Rate	1280 x 1024 pixels Type 160K: 160,000fps (1280 x 8 pixels)	1280 x 1024 pixels Type 200K: 200,000fps (640 x 8 pixels) Type 800K: 800,000fps (640 x 8 pixels) *				
Minimum Exposure	Global electronic shutter 3.9µs (1.01µs at maximum frame rates)					
Dynamic Range (ADC)	12-bit monochrome 36-bit color					
Memory Capacity Options	4GB: 2,180 frames at full resolution 8GB: 4,365 frames at full resolution 16GB: 8,734 frames at full resolution 32GB: 17,472 frames at full resolution					
Memory Partitions	Up to 64 memory segments					
Region of Interest	Selectable in steps of 128 pixels (horizontal) x 8 pixels	kels (vertical) - minimum 640 x 8				
Trigger Inputs	Selectable +/- TTL 5V and switch closure					
Trigger Delay	Programmable on selected input / output triggers: 100ns resolution					
Input / Output	Input: Trigger (TTL/Switch), sync, ready, event, IRIG Output: trigger, sync, ready, rec, exposure					
Trigger Modes	Start, end, center, manual, random					
Time Code Input	IRIG-B					
External Sync	+/- TTL 5Vp-p Variable frequency sync					
Camera Control Interface	High-speed Gigabit Ethernet					
Image Data Display	Frame rate, shutter speed, trigger mode, date/time, status, real time / IRIG time, frame count, resolution					
Saved Image Formats	BMP, TIFF, JPEG, PNG, RAW, RAWW, MRAW, AVI, WMV, FTIF, MOV - Images can be saved with or without image data and in 8-bit, 16-bit or 36-bit depth of sensor where supported					
Supported OS	Microsoft Windows operating system including: XP, Vista, 7, 8, 8.1, 10 (32/64-bit)					

Camera Performance Specifications

* Frame rates above 225,000 fps may be subject to export control regulations in some areas

High-Speed Vigabit Ethernet Interface:

The FASTCALLENSINVIblightamera system is equipped with machingh (Specific Decision Comparison of the system of the

Dedicated I/External Dimensions

A dedicated Bronesa Bediyn for a contact closure hard 20 mm (Hex in 20 m

Camera Body

1.5kg (3.31lbs)

Nikon G-Typer Comparished ens Fitting:

Operational Shock

The FASTCA Operatibly Temperature upped with an object ACO as no diversion of the diversion of the second state of the second

Internal fan cooling (fan-off mode supported) 100G, 10ms, 6-axes



Power


Cooling


Camera Operation Features

Frame Synchronization	Accurate frame synchronization with other cameras and with external and unstable frequencies.		
Memory Partitions	Up to 64 memory segments allow multiple events to be stored in camera memory before downloading, with automatic progression to the next available partition.		
Low Light Mode	Operation at minimum frame rate with separately adjustable shutter time to allow easy camera set-up and focus in ambient lighting.		
IRIG Phase Lock	Enables multiple cameras to be synchronized together with other instrumentation equipment or to a master external time source.		
Internal Time Delay Generator	Allows programmable delays to be set on input and output triggers; 100ns resolution.		
Event Markers	Up to ten user-entered event markers to define specific events within the recorded image sequence .		
Automatic Download The system can be set to automatically download image data to the control PC and, when download is complete to re-arm in readi next trigger with automatically incremented file names.			
Software Binning	Virtual pixel binning (2x2, 4x4 etc.) allows increased light sensitivity with reduced image resolution without changing camera field of view.		

Operation Software Features

Image Calibration	2D image calibration allows the measurement of distance and angle from the image. A calibration grid overlay can be superimposed on the image.
Image Overlay	A stored reference image may be overlaid on the live image to allow accurate camera positioning to achieve the same view as a previous test.
Import of Multiple Image Sequences and Multiple image sequences can be loaded and simultaneously replayed. Timing of image sequences can be adjusted to reference. Time based synchronization allows images captured at different frame rates to be synchronized.	
High Dynamic Range Mode	Making use of the full sensor dynamic range, HDR mode allows enhanced detail in both light and dark areas of an image to be displayed simultaneously.
Motion Detector In order to highlight subtle changes in an image, Motion Detector allows a reference image to be subtracted from a recorded s including propagation of shock waves and surface changes during impact can be visualized using the feature.	
Line Profile A line profile representing grey levels along a line drawn across any region of the image is displayed. In live mode the Line ensure optimum image focus is achieved.	
Histogram	A histogram displaying grey levels within a user-defined image area is displayed. In live mode the Histogram can be used to ensure that optimum exposure levels are set for the scene being recorded.

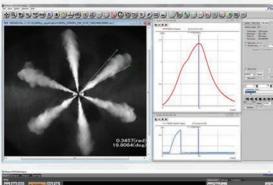

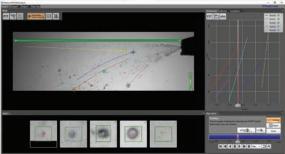


Photron FASTCAM Viewer:

Photron FASTCAM Viewer software (PFV) has been designed to provide an intuitive and feature rich user interface for the control of Photron high-speed cameras, data saving, image replay and simple motion analysis. Advanced operation menus provide access to features for advanced camera operation and image enhancement. Tools are provided to allow image calibration and easy measurement of angles and distances from image data. Also included are a C++ SDK and wrappers for LabView and MATLAB ®.


An optional software plug-in module provides synchronisation between Photron high- speed cameras and data acquired through National Instruments data acquisition systems. Synchronised data captured by the DAQ system provides waveform information which can be viewed alongside high-speed camera images.

Photron FASTCAM Analysis:

PFV software allows image sequences to be exported directly to optional Photron FASTCAM Analysis (PFA) Motion Analysis software. This entry level Motion Analysis software with an on screen 'step by step guide' function launches automatically from Photron FASTCAM Viewer software, and provides automated tracking of up to 5 points using feature or correlation tracking algorithms for the automated analysis of motion within an image sequence.

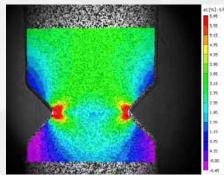


Material Testing & DIC

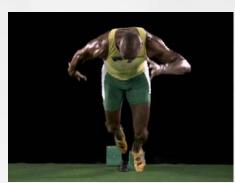

Biomechanics



Welding & Plasma Research

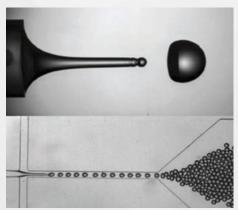

Fluid Dynamics & Microfludics


	Mini UX100						
Resolution	Frame Rate	8GB		16GB		32GB	
(h x v pixels)	Max fps	Frames	Time (sec)**	Frames	Time (sec)**	Frames	Time (sec)**
1280 x 1024	4,000	4,365	1.09	8,734	2.18	17,472	4.37
1280 x 800	6,250	5,587	0.89	11,179	1.79	22,364	3.58
1280 x 720	6,400	6,208	0.97	12,421	1.94	24,849	3.88
1280 x 512	8,192	8,734	1.07	17,468	2.13	34,944	4.27
1280 x 120	40,960	37,248	0.91	74,531	1.82	149,096	3.64
1024 x 1024	4,000	5,458	1.36	10,917	2.73	21,840	5.46
896 x 896	5,120	7,129	1.39	14,259	2.78	28,526	5.57
1280 x 24	204,800	186,242	0.91	372,656	1.82	745,483	3.64

Mini UX50							
Resolution	Frame Rate	8GB		16GB		32GB	
(h x v pixels)	Max fps	Frames	Time (sec)**	Frames	Time (sec)**	Frames	Time (sec)**
1280 x 1024	2,000	4,365	2.18	8,734	4.37	17,472	8.74
1280 x 800	2,500	5,587	2.23	11,179	4.47	22,364	8.95
1280 x 720	3,200	6,208	1.94	12,421	3.88	24,849	7.77
1280 x 512	4,000	8,734	2.18	17,468	4.37	34,944	8.74
1280 x 120	20,480	37,248	1.82	74,531	3.64	149,096	7.28
1024 x 1024	2,000	5,458	2.73	10,917	5.46	21,840	10.92
896 x 896	2,500	7,129	2.85	14,259	5.70	28,526	11.41
1280 x 24	102,400	186,242	1.82	372,656	3.64	745,483	7.28
1280 x 24		186,242		,		,	

opecifications subject to change without notice.

** Recording time is an estimate and may be different depending on recording conditions and settings.

Mechanical & Environmental Specifications

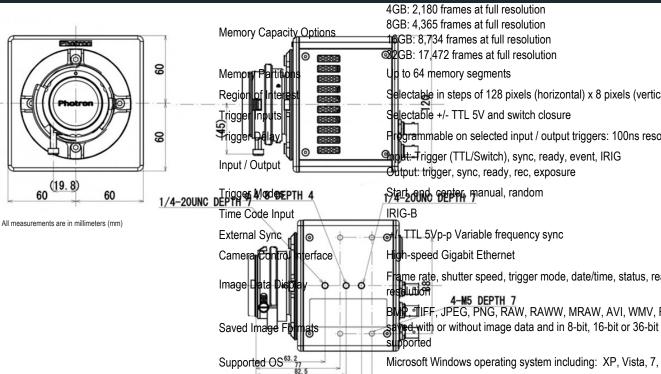

* Frame rates above 225,000 fps may be subject to export control regulations in some areas

Compatibility with Specialist Lens Systems:

A combination of small physical size, low weight and high light sensitivity allows the FASTCAM Mini UX to be coupled to a range of optical systems such as scientific and long distance microscopes, rigid endoscopes or borescopes and image intensifiers for applications ranging from imaging flows in microfluidic devices to combustion diagnostics.

Rugged and Compact Design:

The FASTCAM Mini UX is engineered for use in environments where it may be subject to extreme mechanical shock and vibration. The system has been fully tested for repeated operation at 100G, 10ms, 6-axes.

Small Physical Size:

The small physical size and weight of the Mini camera range allows the use of conventional opto-mechanical hardware for rigid and stable mounting of multiple cameras, and for the location of cameras in space limited locations.

Mechanical and Environmental Specifications

Mechanical	
Lens Mount	F-mount (G-type lens compatible) and C-mount provided
Camera Mountings	4 x 1/4 - 20 UNC (base and top), 4 x M5 (base)
External Dimensions	
Camera Body (excluding protrusions)	120mm (H) x 120mm (W) x 93mm (D) 4.72" (H) x 4.72" (W) x 3.66" (D)
Weight	
Camera Body	1.5kg (3.31lbs)
Environmental	
Operating Temperature	0 to 40C, 32° to 104°F
Storage Temperature	-20 to 60C, -4° to 140°F
Humidity	85% or less (non condensing)
Cooling	Internal fan cooling (fan-off mode supported)
Operational Shock	100G, 10ms, 6-axes
Power	
AC Power (with supplied adapter)	100 to 240V, 50 to 60Hz
DC Power	22 to 32V, 40VA

Specifications subject to change without notice.

PHOTRON USA, INC. 9520 Padgett Street, Suite 110 San Diego, CA 92126 USA

Tel: 858.684.3555 or 800.585.2129 Fax: 858.684.3558 Email: image@photron.com www.photron.com PHOTRON EUROPE LIMITED The Barn, Bottom Road West Wycombe, Bucks. HP14 4BS United Kingdom

Tel: +44 (0) 1494 481011 Fax: +44 (0) 1494 487011 Email: image@photron.com www.photron.com PHOTRON (Shanghai) Room 20C, Zhao-Feng World Trade Building No. 369, JiangSu Road ChangNing District Shanghai, 200050 China Tel: +86 (0) 21-5268-3700 Email: info@photron.cn.com

PHOTRON LIMITED

Kanda Jinbo-cho 1-105 Chiyoda-ku, Tokyo 101-0051 Japan

Tel: +81 (0) 3 3518-6271 Fax: +81 (0) 3 3518-6279 Email: image@photron.co.jp www.photron.co.jp